

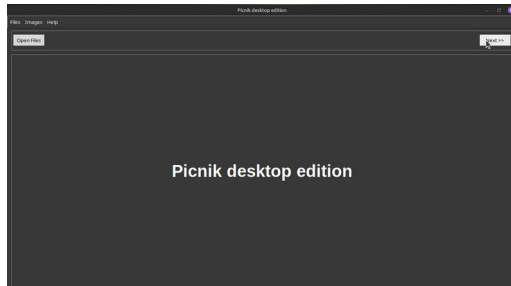
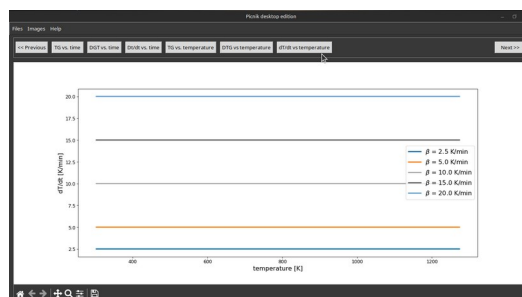
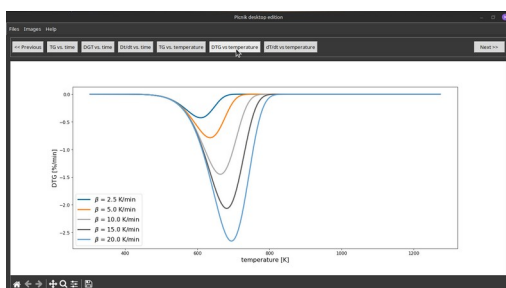
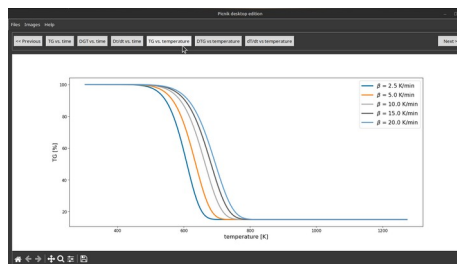
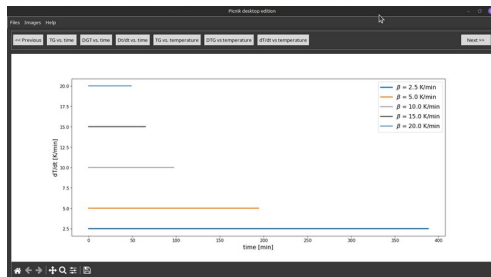


Tutorial picnic Desktop

On the main screen, press the "Open files" button so that the application can read the files, extract the information and generate the first thermogravimetric images that will help the user to enter the correct data in the following steps, then you can go to the next stage by pressing the "next" button and the next screen will be displayed.

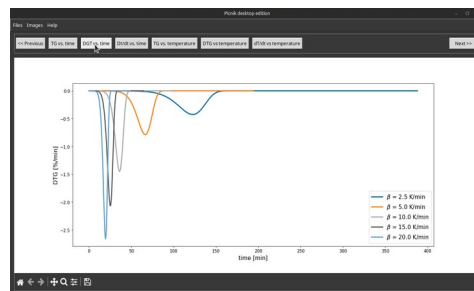


This screen displays the previous and next buttons to return to the previous screen or go to the next, once the user has reviewed the graphs they deem necessary, the user must press the corresponding option to view each graph individually. For example, the "time(min) vs TG(%)" button displays the time(min) vs TG(%) graph.

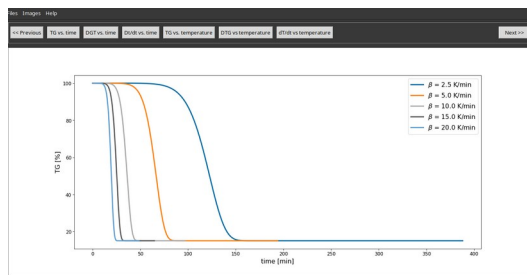

Press the "thermogram" button to view the auxiliary image, then press "enter temperatures" to enter the initial and final temperatures.
The first column is for the initial temperatures and the second for the final temperatures. Once finished, press the "save values" button.





Proceed to save and go to the next step, press next to go to the next screen.


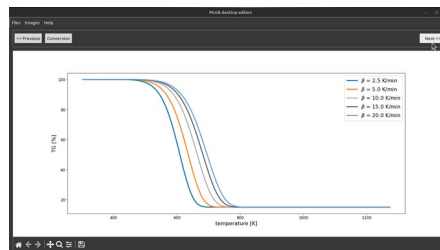






To view the conversion chart, press the "conversion" button; press the "next" button to go to the next screen.

Press the "Isoconversionstep" button to choose the $\Delta\alpha$ you consider appropriate, then press the "next" button to go to the next screen.


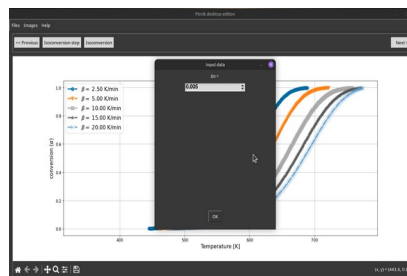
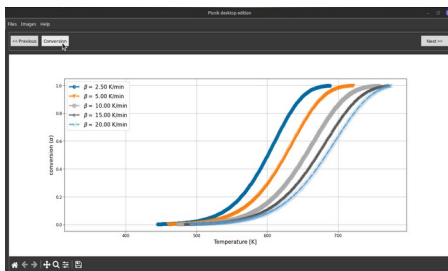


Press the "iso conversion" button to have preview options for time, temperature, and reaction speed, as well as the option to save them. Press next to go to the next screen.



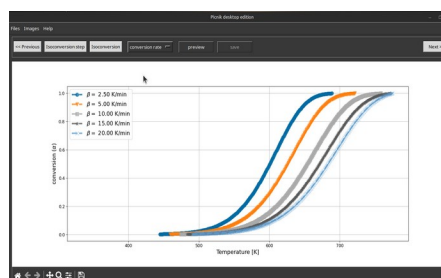
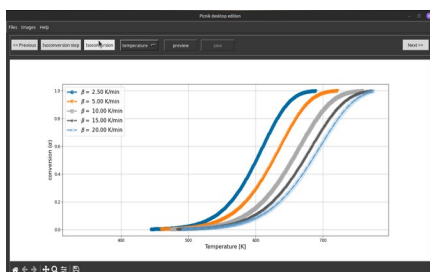


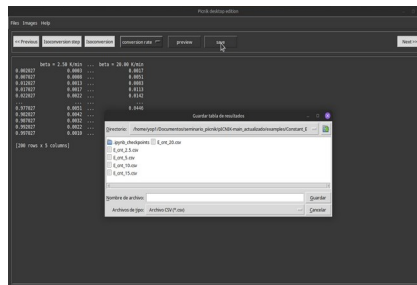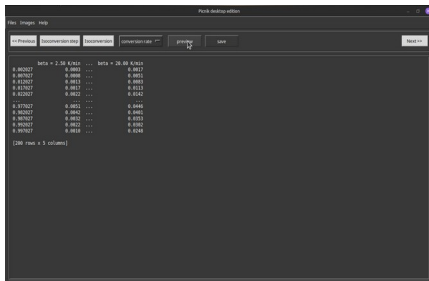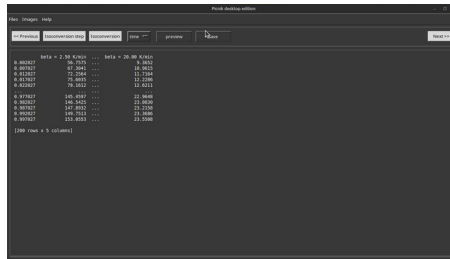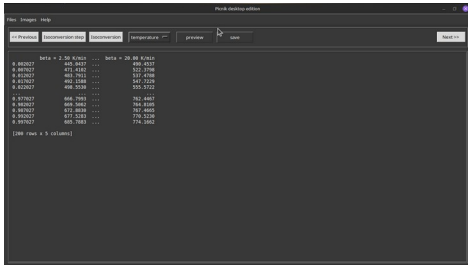


Press the "activation power options" button. A screen will appear with different isoconversion methods. Press the ones you are interested in, then press the "ok" button, and the graphs of the methods you have selected will appear, Press next to go to the next screen.

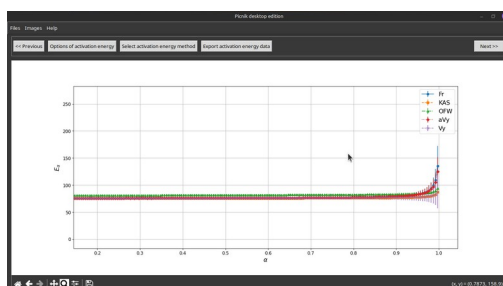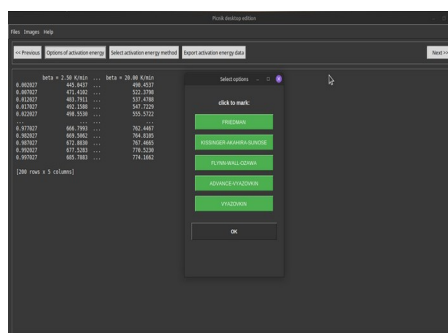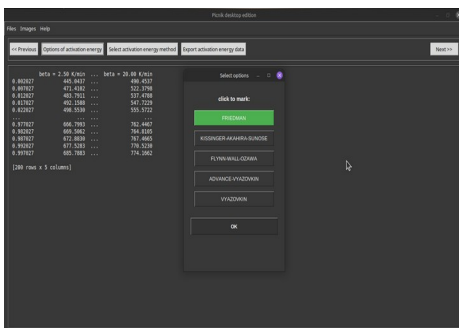


Press the "Choose a method" button and select the method you prefer. You can save your information if you wish. Press "Next" to proceed to the next screen, choose the compensation option you want to view, and you will have the option to save the data, Press next to go to the next screen.

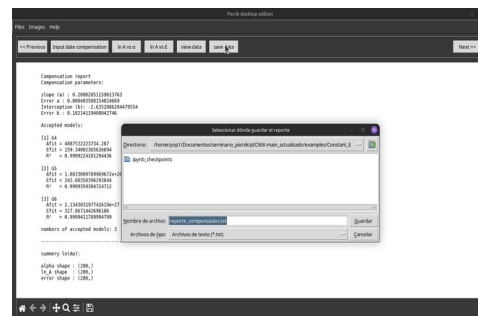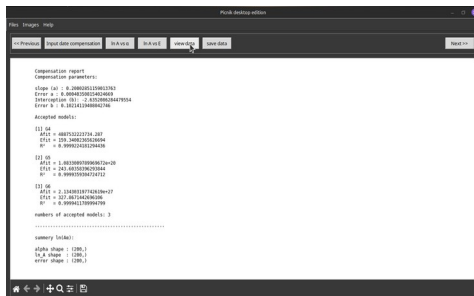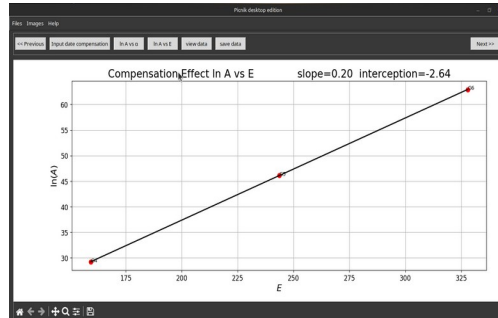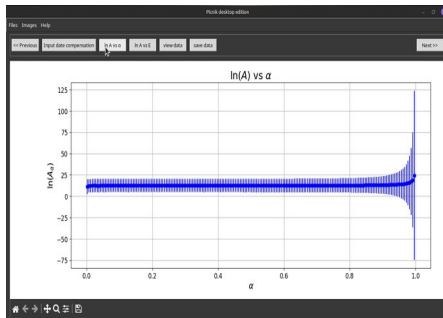

Press the "Reconstruction" button to view the reconstruction graph. Press "Next" to go to the next screen.

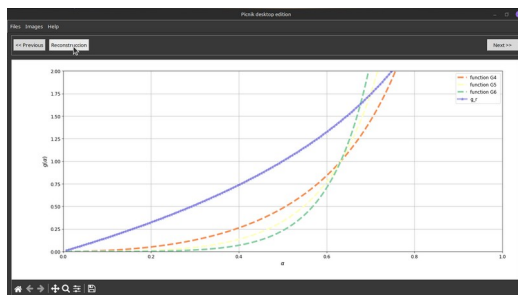

You can export data kinetic, Press next to go to the next screen.

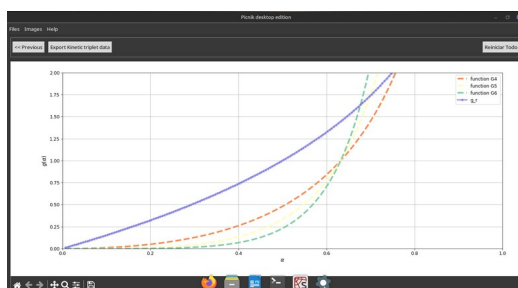

To exit the application, goto menu, select file, and then select exit.

